

SAVE TIME

NO
SMELL

NO
SHRINKAGE

NO
ASH

LC-Resin™



EASY & FAST, LIGHT CURED PATTERN RESIN
(Light Cured JIG GEL)



 **MEDICLUS CO., LTD.**

Head Office & R & D Center :

No.1209, 1210, 1211. 134, Gongdan-ro, Heungdeok-gu, Cheongju-si, Chungcheongbuk-do 28576, Korea
Tel. +82-43-211-2877

Sales Office : No. 735 Brownstone. 464, Cheongpa-ro, Jung-gu, Seoul 04510, Korea / Tel. +82-10-6416-2875

MEDICLUS
MEDI CLUS
MEDICLUS CO., LTD.

LC-Resin™

EASY & FAST, LIGHT CURED PATTERN RESIN

Advantages

- No smell
- No shrinkage
- No ash left in casting process
- Easy to work with optimal viscosity
- Fast working time with light curing type
- Flexural strength : 50Mpa
- Elastic modulus : 1200Mpa
- Curing depth : 5mm

Usage

- Build a jig between implant impression copings
- Build a pattern for bridge, inlay and onlay
- Build a pattern for metal framework of denture
- Build a pattern for clasp of partial denture
- Build a outer shell of telescope

Application



Implant jig






Gold crown



Cap crown

Application process

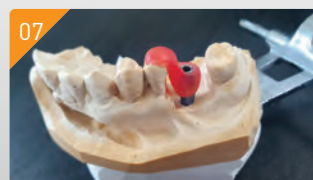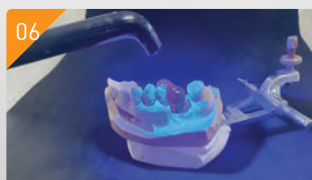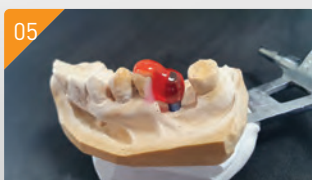

Contents

- 602-1 LC-RESIN RED 14g x 1EA / Disposable tips
- 602-2 LC-RESIN VIOLET 14g x 1EA / Disposable tips

The-Block Out Resin

- 601 The-Block Out Resin 1.2ml x 5syringes
- Disposable tips

